

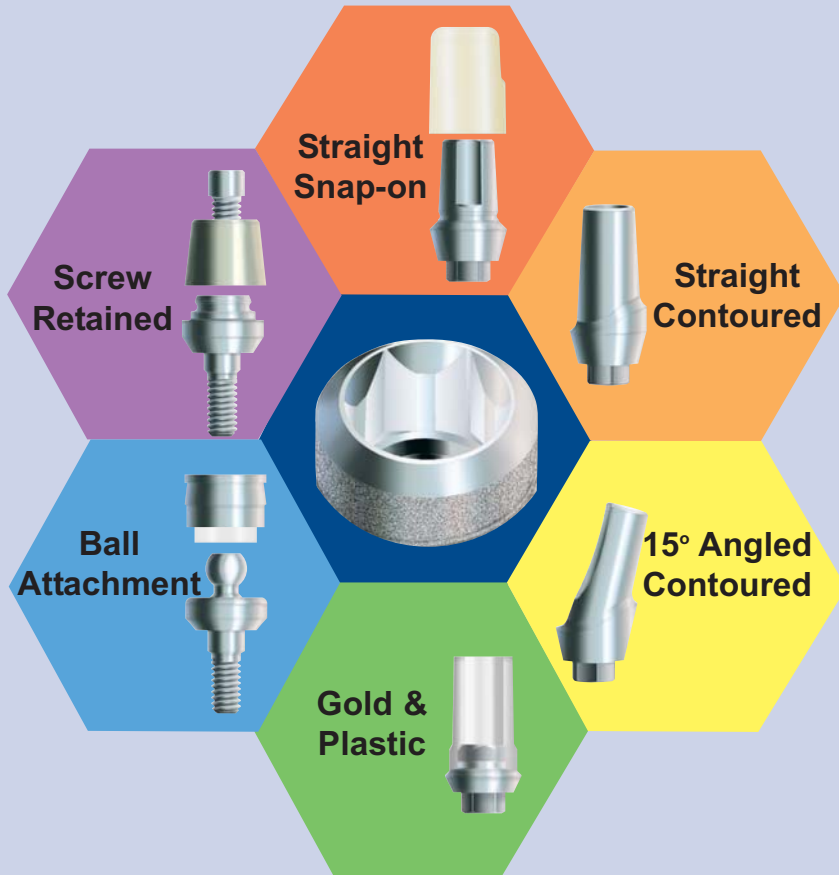


ScrewPlant™ Implant - The Next Generation

From the Developer of the *Screw-Vent® Implant

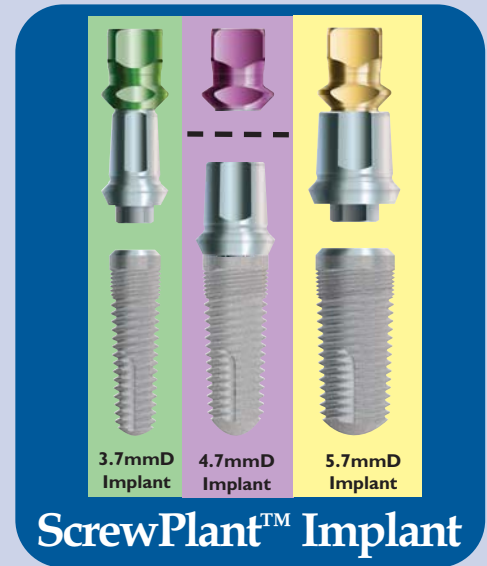

Shorten
Fixture-mount
for Final
Abutment
Snap-on
Comfort Cap
Included

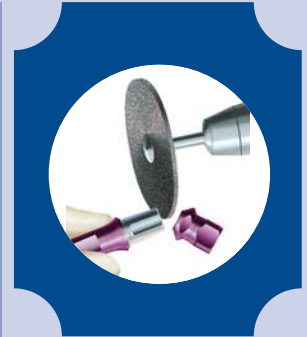

Implant Direct also offers Legacy Abutments Compatible with Zimmer Dental's Screw-Vent and BioHorizon's Prodigy Implants
 Contoured Abutments: *Zimmer Dental \$145-\$165; Implant Direct \$85
 * Screw-Vent® and Core-Vent® are registered trademarks of Zimmer Dental
 ** Prodigy Implant is a registered trademark of BioHorizons Inc.

ScrewPlant™ Implant
 25 Years of Innovation

1982 1986 1999


*Core-Vent® > *Screw-Vent® > *Tapered SV®

Compare ScrewPlant™ to *Screw-Vent®

- Save 70% - Healing Collar & Abutment Included
- Platform: External Bevel for Stronger Interface
- Surgically Compatible with *Screw-Vent® Drills
- SBM (HA Blasted) surface to top of implant
- Mini-threads for reduced stress near crest
- Eliminate Vent and lengthen vertical cutting groove
- Even Taper over entire body for bone expansion

SURGICAL COMPATIBILITY WITH



<p>ScrewPlant™ [70% Savings] U.S. List Price = \$150 Includes 2mm Healing Collar & Abutment</p>	<p>Zimmer Screw-Vent® U.S. List Price Implant = \$513 for Implant(\$325), HC(\$43), Abutment (\$145)</p>
--	---

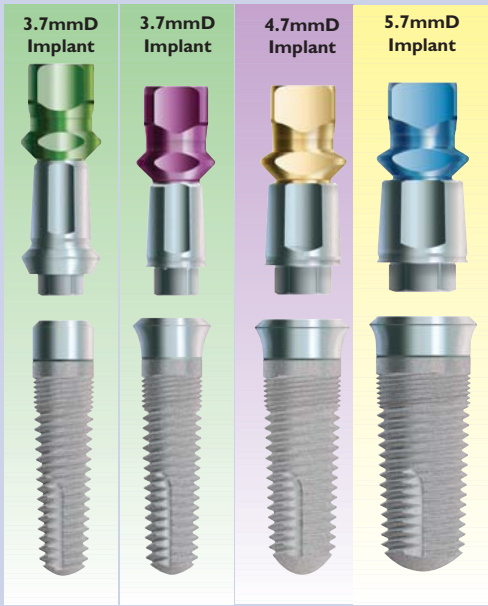



ScrewPlus™ 1-Stage Implant

With All-in-one Packaging
and SBM Surface
with 15 years of Proven
Clinical Success


**ScrewPlant and ScrewPlus
Accept Same Abutments**

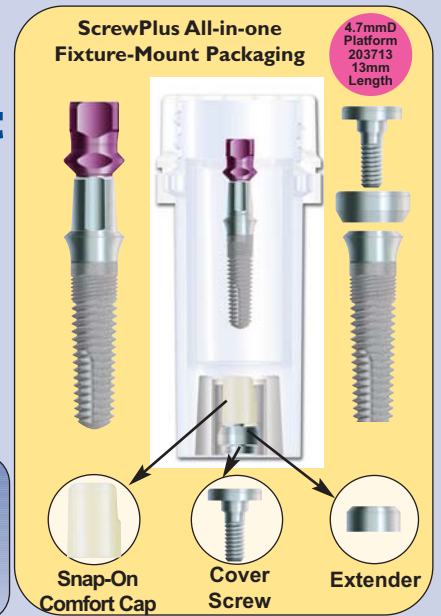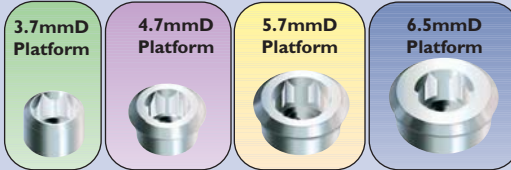




IMPLANT DIRECT'S SBM SURFACE - A DECADE OF CLINICAL SUCCESS

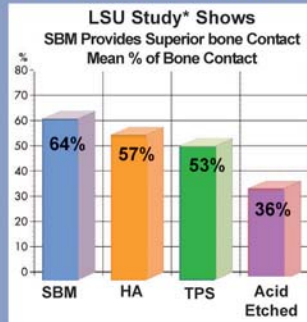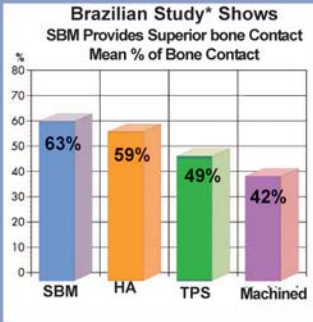

SBM® Soluble Blast Media
A medium-rough surface created
by blasting with soluble HA,
eliminating the need for a
secondary acid etch procedure that
further rounds the tips of the threads

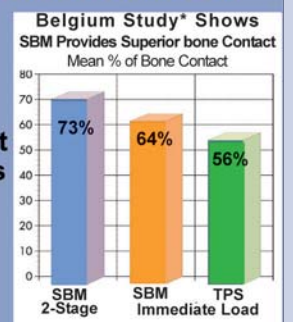

First introduced with the *Screw-
Vent implant in 1992 and used
today on Zimmer, BioHorizons,
Lifecore & Implant Direct

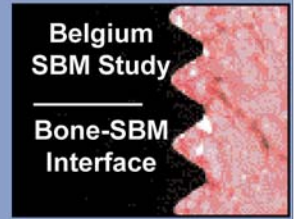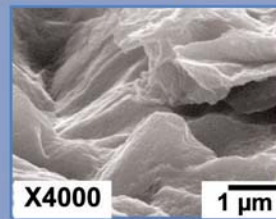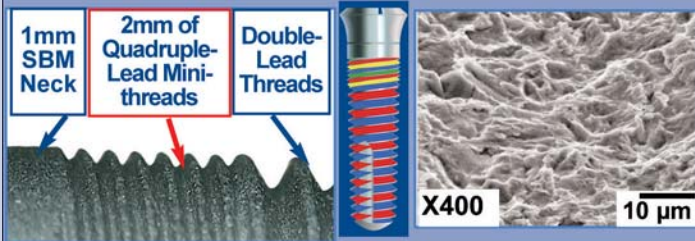

ScrewPlus™ Implant

Compare ScrewPlus™ to *ITI®

- Save 74% based on U.S. List Prices
- Packaged on color-coded Fixture-mount/Transfer
- Fixture Mount Modified to Final Abutment
- Implant vial includes 2mm Healing Collar
- 2mm deep Internal Hex vs 0.5mm deep Octagon
- Surface: Blasted with Soluble HA (SBM™)
- Mini-threads for reduced stress near crest
- Double-lead threads for faster insertion
- Even Taper for Soft Bone Compression
- Vertical self-tapping, cutting groove + Apical Threads

PRICE AND DESIGN ADVANTAGES COMPARED TO

Implant Direct ScrewPlus™	Straumann Tapered ITI®	Compare U.S. List Prices	ITI Tapered Implant	Screw Plus Implant
		Implant	\$345.00	\$150
		Cover Screw	\$34.50	Included
		Healing Cap	\$48.00	2mm Included
		Straight Abutment	\$131.00	Included
		Snap-on Cap	\$7.50	Included
		Snap-on Transfer	\$7.50	Included
		Total US List Price	\$573.50	\$150
		SAVINGS	3.8X	74%

*ITI® is a Registered Trademark of Straumann