

Gerald A. Niznick, DMD, MSD  
 President, Implant Direct LLC  
 Inventor of 30 Dental Implant Related  
 Patents including the Internal  
 Anti-rotational Connection  
 U.S. Patent #4,960,381




**Implant Direct's Mission Statement:**  
 Provide experienced implant dentists with  
 high quality, innovative implants and  
 abutments through internet sales at  
 factory-direct prices.

## Spectra-System™ Application Specific Implants

The Spectra-System™ product line is the fulfillment of Dr. Niznick's vision to create an implant system with a universal body design and a variety of abutment platforms, providing a spectrum of clinical applications. All-in-one packaging is a uniform feature of the five Spectra-System implants, providing the abutment, transfer, healing collar and/or comfort cap for one low price. Each tapered implant body features 2mm of "mini-threads" near the top of the implant which have been known to reduce stress concentration near the crest of the ridge\*, and all but the ScrewRedirect 1-piece implant provide double lead threads for faster insertion. The ScrewPlant™ 2-stage implant and ScrewPlus™ 1-stage implant both have the same external bevel/internal hex platform. These 2-piece implants are each provided on a color-coded fixture-mount that is used as a transfer and can be shortened to create a straight abutment with a 1mm collar for cemented restorations. These implants are also packaged with a cover screw, 2mm healing collar and a snap-on comfort cap for the modified straight abutment. The ScrewDirect™ 1-piece implant has a 1.5mm collar and the same straight abutment head (beveled shoulder/retentive groove) as that created by shortening the fixture mounts. It is packaged with a snap-on comfort cap and color coded snap-on transfer that matches the color of an aluminum abutment analog. The Redirect™ 1-piece implant provides a 12 degree angled head with a gingival margin that is contoured from 0.75mm on the labial to 2.5mm on the lingual. It is indicated for cemented crowns in the maxillary anterior region. The Indirect™ 1-piece implant has a screw-receiving platform for bars or fixed-detachable restorations. It has a 2.5mm collar height and is packaged with a 2mm extender, providing up to 5mm in total supra-crestal height. Also included is a snap-on comfort cap and transfer that engage a retentive groove above the 5mmD wide platform. US List Price for all 5 implants is \$150 (2007-8).

### One-Piece Spectra-System Implants

ScrewDirect™ 1-Piece 1.5mm Neck	ScrewDirect™ 1-Piece 1.5mm Neck	ScrewIndirect™ 1-Piece 2.5mm Neck	ScrewRedirect™ 1-Piece 1.0mm Neck
---------------------------------------	---------------------------------------	-----------------------------------------	-----------------------------------------

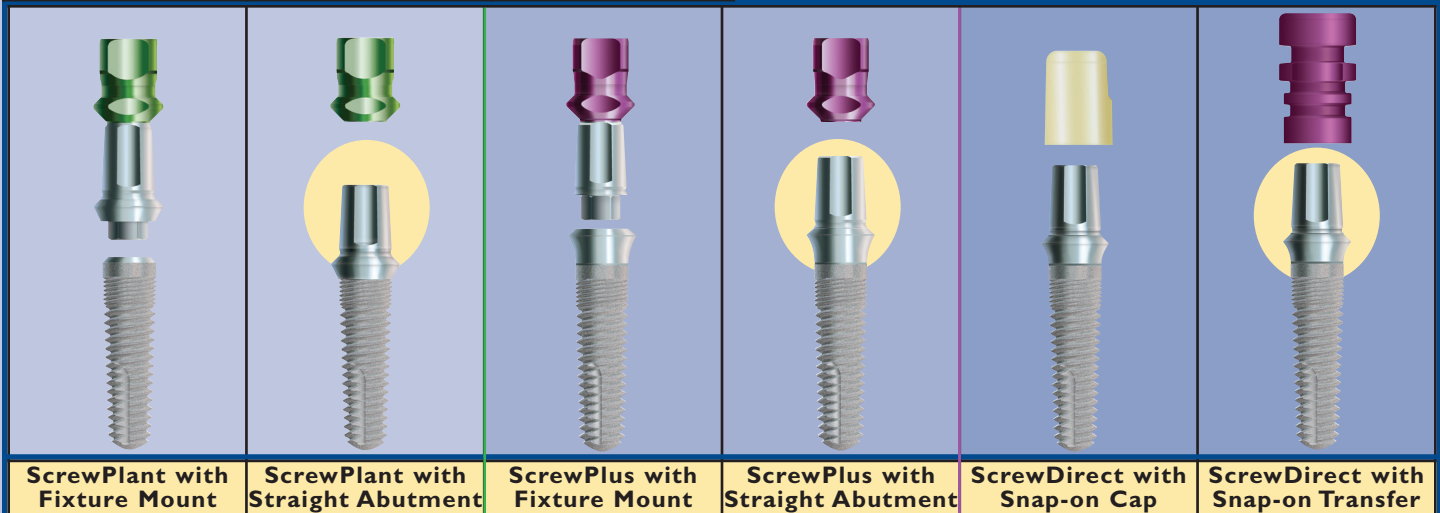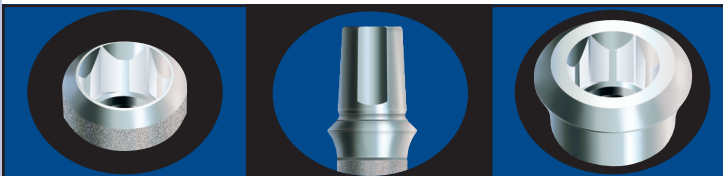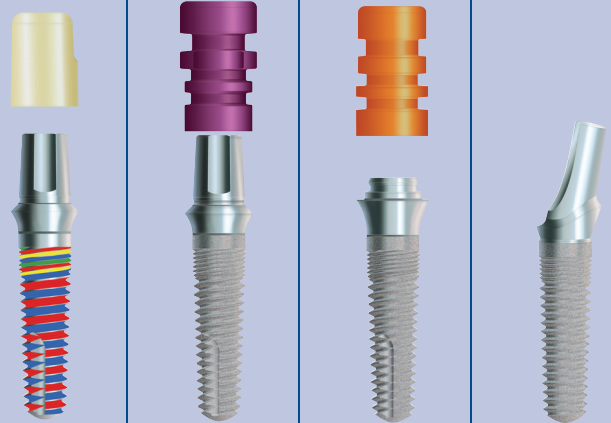

ScrewPlant with  
Fixture Mount

ScrewPlant with  
Straight Abutment

ScrewPlus with  
Fixture Mount

ScrewPlus with  
Straight Abutment

ScrewDirect with  
Snap-on Cap

ScrewDirect with  
Snap-on Transfer



27030 Malibu Hills Road  
 Calabasas Hills, CA 91301  
 Phone: 818.444.3333 Fax: 818.444.3400  
 www.implantdirect.com

\* Stig Hansson. The Implant Neck: Smooth or Provided with Retention Elements. A Biomechanical Approach.  
 Clin Oral Impl Res 1999;10:394-405