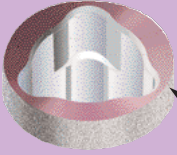
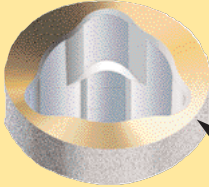
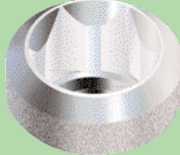
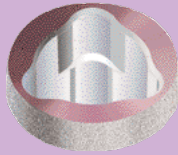
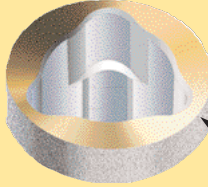

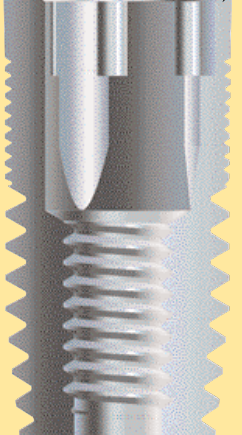
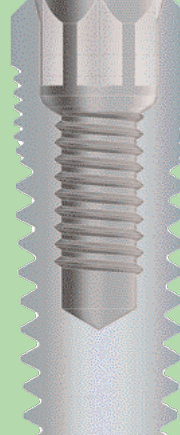
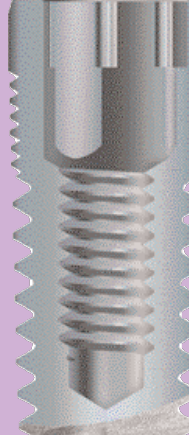
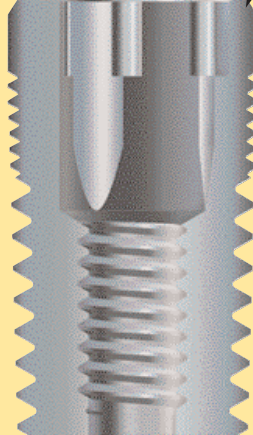









3.5mmD Platform	4.3mmD Platform	3.7mmD Platform	3.5mmD Platform	4.3mmD Platform
				
Wall Thickness 0.009"	Wall Thickness 0.012"	Wall Thickness 0.018"	Wall Thickness 0.013"	Wall Thickness 0.020"
				
Replace/RePlant 3.5D	Replace/RePlant 4.3D	ScrewPlant 3.7D	RePlus 3.7D	RePlus 4.7D

## The RePlus Implant uses Implant Direct's Drills and Surgical Protocols

 3.5mmD Platform  4.3mmD Platform  5.0mmD Platform  3.7mmD  4.7mmD  5.7mmD 	1	<b>50% stronger Small Diameter Implant</b>
	2	<b>Micro-threads for reduced stress near crest</b>
	3	<b>Double-lead threads for 2 times faster insertion</b>
	4	<b>Self-tapping vertical cutting groove 1/2 the length</b>
	5	<b>Soft Bone Protocol: Evenly tapered for Expansion</b>
	6	<b>Hard Bone Protocol: No need for bone tap</b>
	7	<b>Threaded to the apex for increased stability</b>
	8	<b>Straight, step-drills: faster surgery with less drills</b>
	9	<b>Five Implant Length Options, including 11.5mm</b>
	10	<b>Two Surface Options: HA Coated or Blasted</b>
	11	<b>Precision Manufacturing: &lt;math&gt;&lt;0.5^\circ&lt;/math&gt; vs &lt;math&gt;&gt;1.5^\circ&lt;/math&gt; Rotation</b>
	12	<b>All-in-1 Pkg. for 70% *Savings: \$200 vs \$641</b>

The RePlus Implant is evenly tapered over its entire length, allowing insertion into an undersized socket in soft bone. Research has shown this not only increases initial stability essential for immediate loading, it also results in increased bone attachment and removal torque following osseointegration. In hard bone, Implant Direct's final step-drills allow insertion without the need for a bone tap. RePlus Implants are available with SBM medium rough surface over their entire length and with HA coating starting 3mm from the top of the implant. The Trilobe, color coded internal connection accepts Nobel's Replace and Implant Direct's RePlant abutments. Each RePlus vial includes a cover screw, titanium Fixture-mount/transfer that can be shortened for use as the final "snappy" abutment. A snap-on comfort cap is also provided. The 3.7mmD RePlus implant provides 50% greater strength than NobelRePlace 3.5mmD Implant. The transition from 3.7, 4.7 and 5.7mmD bodies to the 3.5, 4.3 and 5.0mmD Trilobes create platform switching.

\*Includes Implant, Abutment, Transfer & CS

