



New Legacy™ Angled Overdenture Abutments Expand Treatment Options

New Legacy™ Angled Overdenture Abutments provide unmatched versatility at approximately half the cost of comparable abutments for greater overdenture treatment plan flexibility and case acceptance.

Compatible with nine implants across the Legacy™ System, these 15° and 30° angled abutments allow you to treatment plan with an assortment of body types and diameters that utilize a single surgical kit. You can even choose among the numerous internal-hex implants, ranging from 3.3mm to 6.0mm diameters, that these hex-engaging abutments are compatible with industry-wide.

In addition to screw-retained prostheses, the angled abutment base may be used with the your preferred overdenture attachment system when combined with a ball top or a GPS™ abutment top featuring a LOCATOR®-compatible profile (sold separately).

With the “simply smarter” design, you can easily seat the abutment screw without comprising the integrity of the oral cavity seal – unlike some current market offerings.

The 15° and 30° angled abutments are available for both 3.5mm and 4.5mm platforms. Our All-in-One packaging includes an abutment screw, transfer, screw-receiving top and comfort cap along with the abutment base.



Expand your overdenture treatment options today!

[Watch Tutorial](#)

[Shop Online Now](#)



All-in-One packaging includes angled abutment base, abutment screw, screw top, transfer and comfort cap.

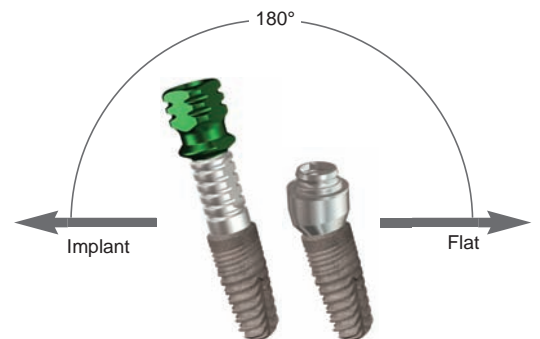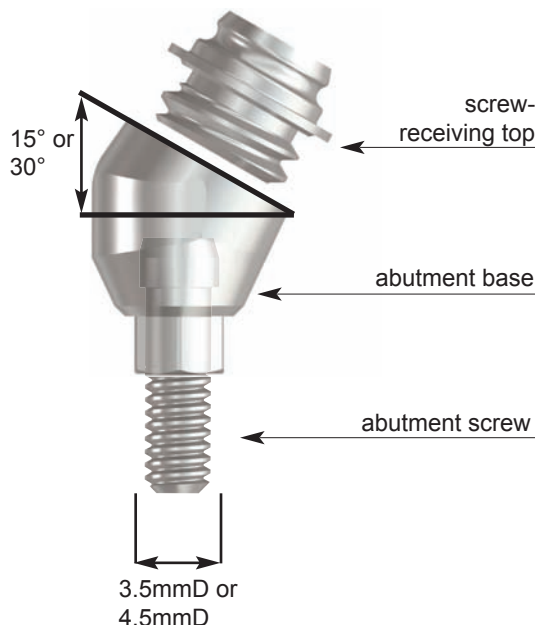

During implant placement, the longitudinal flat of the fixture mount **must** be oriented 180° from implant angulation