



Swish™ System Catalog



SwishPlus™

SwishPlant
with Transfer
& Abutment



Tapered with Micro-Threads
& Double-Lead Threads

SwishPlus
with Straumann®
Compatible Carrier




Straight with Micro-Grooves
& Single-Lead Threads

Designed For Success, Packaged For Value

Customer Service & Technical Support online, over the phone, or in your office

Five Year Warranty²

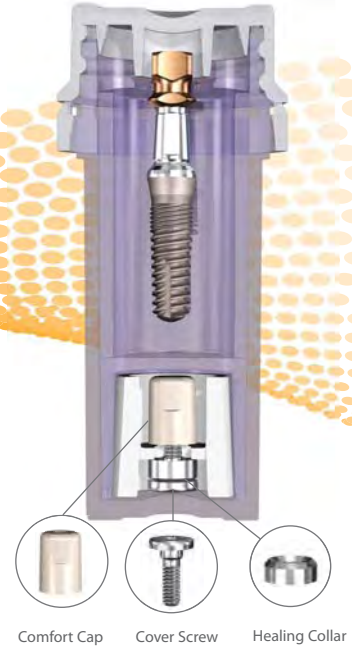

Reality Check
Straumann Customers
Save \$499 with SwishPlant

SwishPlant™ Implant

All-in-1 Packaging includes implant, straight abutment/transfer, cover screw & healing collar
\$200 vs \$699¹ from Straumann®

Swish™ Implant System Advantages:

Internal Octagon (Niznick US Pat. #4,960,381 Licensed to Straumann)
Prosthetic compatibility with Straumann's Tissue-Level implants (Standard, Standard Plus and Tapered Effect)
Internal square broached in flats of octagon for fixture-mount engagement (Pat. Pend.)

Two Implant Designs & Packaging Options

Allows for selection based on price, packaging or thread design

SwishPlus³: \$150 includes cover screw, healing collar & carrier compatible with Straumann tool & standard ratchet

- 1 2mm neck half textured for tissue-level or esthetic placement
- 2 Micro-grooves for improved tissue attachment & reduced stress
- 3 Progressively-deeper single-lead threads for increased surface area & complete surgical compatibility
- 4 Self-tapping groove for increased initial stability
- 5 Tapered apical end to initiate self-tapping

SwishPlant: \$200 includes cover screw, healing collar, comfort cap & straight abutment/transfer


- 1 2mm neck lightly textured for bone- or tissue-level placement
- 2 Micro-threads for improved tissue attachment, reduced stress & increased stability
- 3 Progressively-deeper double-lead threads for increased surface area & faster insertion
- 4 Self-tapping groove for increased initial stability
- 5 Tapered body for soft bone expansion

SwishPlus

Straight with Micro-Grooves & Straumann Threads

SwishPlant

Tapered with Micro-Threads & Double-Lead Threads

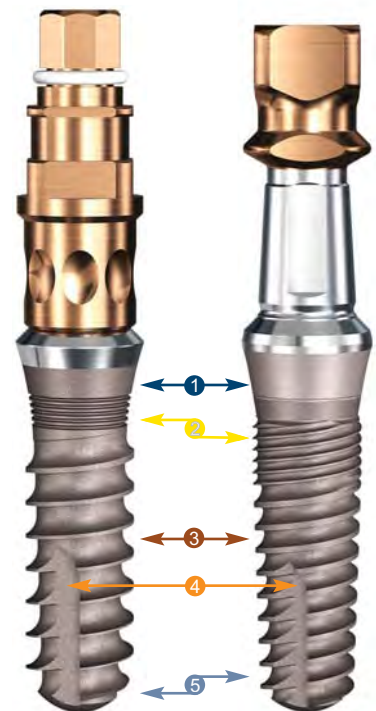

	Solid Abutment Kit	Straight Narrow Abutment	Straight Snap-On	Straight Contoured	15° Angled Contoured	Gold/Plastic	Zirconia/Ti Abutment	Angled Zirconia/Ti Abutment	Plastic Temporary Abutment	Ball Attachment	GPS™ Attachment	Screw Receiving
Our price	\$100	\$85	\$85	\$85	\$85	\$100	\$120	\$120	\$35	\$107	\$100	\$85
Straumann	\$238	\$181	\$181	\$182	\$203	\$157	\$339 ⁴	\$339 ⁴	\$61	\$192	\$193	\$169

Intro Offer: Make the switch & receive three FREE implants?

www.implantdirect.com | 888-649-6425

Price comparisons based upon US list prices as of April 2011. All trademarks are property of their respective companies. ¹US list price for SLActive Tapered Effect implant, closure screw, healing abutment, solid abutment, transfer and comfort cap. ²Terms and conditions apply. ³System Variation: 3.3mmD & 5.7mmD SwishPlus are tapered; 3.3mmD also has SwishPlant threads. ⁴Comparison based upon Straumann's CARES Custom Zirconia Abutment. Price varies depending upon the laboratory.



7K-09-11-0033 Printed in the U.S.A.

Copyright © 2011 Implant Direct Sybron International. All Rights Reserved