

3.2mmD Straight /Tapered Implant

3.0/3.2



Soft Bone Final Dense Bone Final Optional Bone Contour



SD2.3
2.3/2.0mmD



SD2.8
2.8/2.3mmD



CBD3.2

3.7mmD Straight /Tapered Implant

3.5/3.7



SD2.3
2.3/2.0mmD

Soft Bone Final Dense Bone Final Optional Bone Contour



SD2.8
2.8/2.3mmD



SD3.4
3.4/2.8mmD



CBD3.7

4.2mmD Straight /Tapered Implant

3.5/4.2



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD

Soft Bone Final Dense Bone Final Optional Bone Contour



SD3.4
3.4/2.8mmD



SD3.8
3.8/3.4mmD



CBD4.2

4.7mmD Straight /Tapered Implant

4.5/4.7



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD



SD3.4
3.4/2.8mmD

Soft Bone Final Dense Bone Final Optional Bone Contour



SD3.8
3.8/3.4mmD



SD4.4
4.4/3.8mmD



CBD4.7

5.2mmD Straight /Tapered Implant

4.5/5.2



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD



SD3.4
3.4/2.8mmD

Optional



SD3.8
3.8/3.4mmD

Soft Bone Final Dense Bone Final Optional Bone Contour



SD4.4
4.4/3.8mmD



SD4.8
4.8/4.4mmD



CBD5.2

5.7mmD & 7.0mmD Straight /Tapered Implant

5.7/5.7

5.7/7.0



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD

Optional




SD3.4
3.4/2.8mmD

Optional



SD3.8
3.8/3.4mmD

Optional

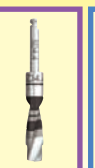

SD4.4
4.4/3.8mmD

Optional



SD4.8
4.8/4.4mmD

Soft Bone Final



D5.1
5.1mmD

Dense Bone Final



SD5.4
5.4/4.8mmD

Optional Bone Contour for 5.7mmD



CBD5.7

Optional Bone Contour for 7.0mmD



CBD7.0