

Dr. Travis Perkins  
770 West Sahara # 114  
Las Vegas, NV 89117

Attn: Dr. Travis Perkins

February 03, 2010  
P.O. #: 11819  
AMENDED REPORT  
Original Issue Date:  
01-22-10  
❖ Amendment Summary

## PACKAGING TEST REPORT

**Sample Information:** Dye Immersion Test on Nobel Biocare Implants  
1. Ref. 32225 Lot # 400832 (control)  
2. Ref. 32225 Lot # 400832  
3. Ref. 32214 Lot # 401060  
4. Ref. 32212 Lot # 388864  
5. Ref. 32218 Lot # 678113  
6. Ref. 32218 Lot # 674453  
7. Ref. 32213 Lot # 410954  
8. Ref. 32216 Lot # 400709

**Date Received:** January 15, 2010  
**Date in Test:** January 20, 2010  
**Date Completed:** January 20, 2010

**Test Information:** Test Code: 38091  
Procedure #: PK203SOP, PK203IMD.201  
Gauge of Needle Used: N/A  
Pressure (mmHg): 580  
Dwell Time: 1 minute  
Visible Light (Lx): >1000 Lux  
Lux Meter ID / Calibration Due: MTR0004 / 09-10  
Dye Lot # / Expiration Date: 2421 / 04-07-10  
Time ID / Calibration Due: TMR0032 / 04-10

### Sample 1

Sample ID	Results
1 - Control	FAIL
2	FAIL
3	FAIL
4	FAIL
5	FAIL

(continued)



Report Number  
**827731.A01.R01**  
Page 2 of 2

Dr. Travis Perkins  
770 West Sahara # 114  
Las Vegas, NV 89117

Attn: Dr. Travis Perkins

February 03, 2010  
P.O. #: 11819  
AMENDED REPORT  
Original Issue Date:  
01-22-10  
❖ Amendment Summary

## PACKAGING TEST REPORT

### Sample 1

Sample ID	Results
6	FAIL
7	FAIL
8	FAIL


❖ A01 – Dated 01-26-10: Sponsor modified – Company name.

*Testing conducted in accordance with current Good Manufacturing Practices.*





Before dye immersion testing

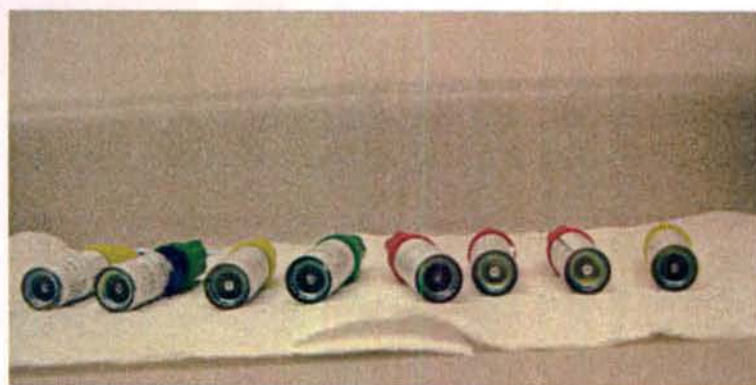

After testing (control is 2<sup>nd</sup> from left)



After testing (control is 2<sup>nd</sup> from left)